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Product Name


Pressure Therapy System Beauty Line Manual.

Use

The products of the Beauty Line series are devices for lymphatic drainage; a treatment against water retention. The pressure starts from the distal section and moves in the proximal section of the leg cuff. In this way, the flow of body fluids is stimulated, facilitating their elimination. The pressure therapy is a medical device, do not use without a physician opinion in case of health problems.

Pump Features	
Code	07.01.HF3006
Dimensions	24 x 18 x 11 cm
Display	Manual
Pressure Range	20-250 mm Hg
Cycle Time	From 10 to 30 min.
Power	220/240 V - 50 Hz
Consumption	30 Watt
Noise Level	≤ 35 Dba
Massagge Mode	A • From Distal to Proximal

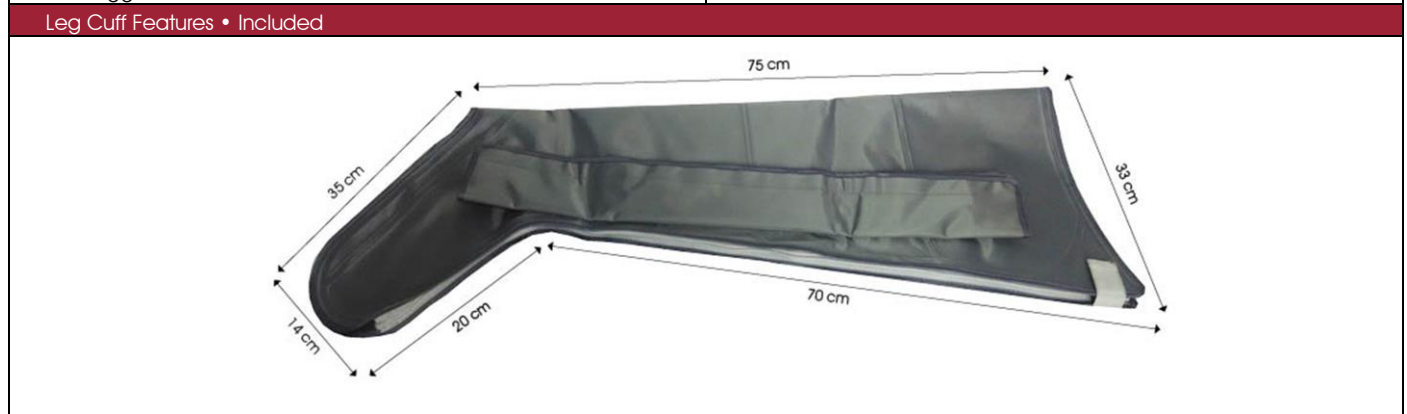Material	Nylon • PU
----------	------------

Overlapping model: the cells that make up the leg cuff are overlapped. This feature allows for a more gradual and effective compression

Warnings
Pressure therapy is contraindicated for people suffering from the following pathologies: cardiopathy, suspected or ascertained thrombosis, arteriosclerosis or ischemia, dermatitis, gangrene, local infections
Check the correct assembly of the device with particular attention to the connection of the tubes between the pump and the leg cuff

Maintenance and Cleaning
Clean the leg cuff with a neutral soap; the drying must not take place in direct sunlight
Do not immerse the tubes or the pump in water
Check the connection tubes are not obstructed by dust. In this case use a soft swab to clean them

Storage and Disposal
Single-user Product
Keep away from heat, sparks, open flame or other sources of ignition
Dispose of the device in compliance with current regulations on environmental protection and separate waste collection

Package and Nr. Items	
Box	1 Pump, 1 Leg Cuff Set • 2 Pcs

Warranty	
European	24 months for manufacturing defects with the exception of consumable parts

Product References	
Manufacturer	Demarta Virginio S.r.l.
Directive MDR 2017/745	Class I
CND Code	Z120607
RDM Code	1745140/R